

Diagnostic X-Ray Shielding

Radiographic/Fluoroscopic Rooms Multi-Slice CT Rooms Using NCRP 147 Methodology

Melissa C. Martin, M.S., FAAPM, FACR
Therapy Physics Inc., Torrance, CA

AAPM Annual Meeting, Orlando, FL
Refresher Course
Thursday, July 26, 2007 7:30 am

Acknowledgement

**Radiographic /Fluoroscopic Room
Slides Courtesy of:**

Ben Archer, Ph.D, FACR, FAAPM
Baylor College of Medicine, Houston, TX

Required Information for Shielding Designs

- Architectural drawings of equipment layout in room
- Architectural drawings of surrounding areas indicating usage of these areas - offices, restrooms, corridor, exterior, etc.
- Elevation view of room or construction of floor and ceiling and distance between floors

Nomenclature for Radiation Design Criteria

Required thickness = NT/Pd^2

where:

N = total no. of patients per week

T = Occupancy Factor

P = design goal (mGy/wk)

d = distance to occupied area (m)

Shielding Design Goal (Air Kerma):

Uncontrolled Areas

Annual: $P = 1$ mGy per year
Weekly: $P = 0.02$ mGy per week

Controlled Areas

Annual: $P = 5$ mGy per year
Weekly: $P = 0.1$ mGy per week

New Formalism for Radiation Design Criteria

Required thickness = NT/Pd^2

where:


N = total no. of patients per week

T = Occupancy Factor

P = design goal (mGy/wk)

d = distance to occupied area (m)

Easy to use graphs for R and RF rooms developed by Simpkin are included in Report.

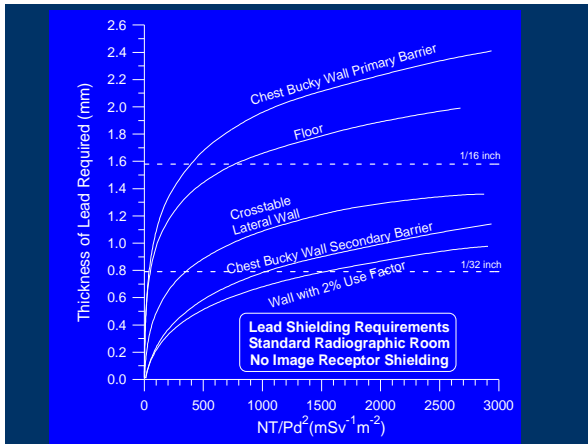

Shielding Design Goal (Air Kerma):

Uncontrolled Areas

Annual: $P = 1$ mGy per year
Weekly: $P = 0.02$ mGy per week

Controlled Areas


Annual: $P = 5$ mGy per year
Weekly: $P = 0.1$ mGy per week

Distance (d)

The **distance in meters** from either the primary or secondary radiation source to the occupied area.

New recommendations in Report 147 for areas above and below source.

Where in the occupied area do you calculate the dose?

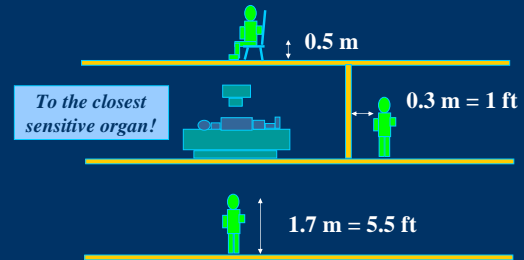

Recommended Occupancy Factors for Uncontrolled Areas:

$T=1$ Clerical offices, labs, fully occupied work areas, kids' play areas, receptionist areas, film reading areas, attended waiting rooms, adjacent x-ray rooms, nurses' stations, x-ray control rooms

$T=1/2$ Rooms used for patient examinations and treatments

$T=1/5$ corridors, patient rooms, employee lounges, staff rest rooms

$T=1/8$ **corridor doors**

Recommended Occupancy Factors for Uncontrolled Areas:

$T=1/20$ public toilets, vending areas, storage rooms, outdoor area with seating, unattended waiting rooms, patient holding areas

$T=1/40$ **minimal occupancy areas**; transient traffic, attics, unattended parking lots, stairways, janitor's closets, unattended elevators

Pre-shielding (x_{pre}) for Radiographic Room Workload Distributions

(Dixon RL, Med Phys 1994)

Grid + cassette: (cross table)
 Equivalent to: 0.3 mm Pb
 or 3 cm concrete

Grid + cassette + table/chest bucky supports: (over table and chest)
 Equivalent to: 0.85 mm Pb
 or 7.2 cm concrete

Equivalency of Shielding Materials

Table 4.8 Page 67

Steel thickness requirement:

8 × Pb thickness requirement

Gypsum wallboard thickness requirement:

3.2 × concrete thickness requirement


Plate Glass thickness requirement:

1.2 × concrete thickness requirement


Light-weight concrete thickness requirement:

1.3 × std-weight concrete thickness requirement

Figure 5.2 Radiographic Room page 75

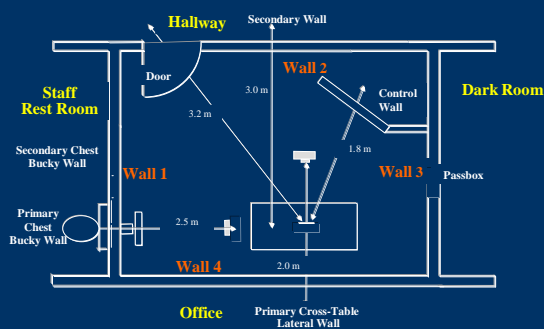

PRIMARY BARRIER Cross-Table Wall in Rad Room

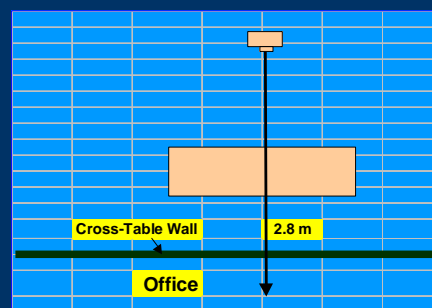

Simplified Graphical Solution Cross-Table Wall in Rad Room

Required thickness $\rightarrow NT/Pd^2$

where:

N = 125 patients/ week

T = 1

P = 0.02 mGy/wk


d = 2.8 m

$$NT/Pd^2 = 797 \text{ mGy}^{-1} \text{ m}^{-2}$$

Simplified Graphical Solution Cross-Table Wall in Rad Room

1. Go to page 54, Fig. 4.5a
(Primary, lead, with no pre-shielding)
2. Look up $NT/Pd^2 = 797$
(Cross-table Wall)

Pb required = 1.03 mm
Specify: 4/64" (1/16"); 4 lb/sqft

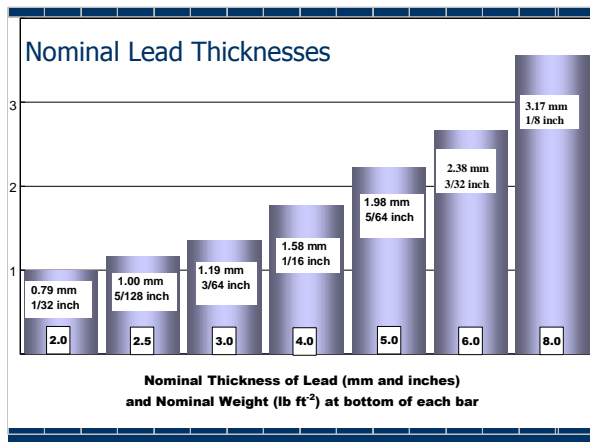

Simplified Graphical Solution Cross-Table Wall in Rad Room

OR


1. Go to page 55, Fig. 4.5b
(Primary, lead, with pre-shielding)
2. Look up $NT/Pd^2 = 797$
(Cross-table Wall)

Pb required = 0.83 mm
Specify: 5/128"; 2.5 lb/sqft (minimum)
Recommended: 1/16"; 4 lb/sq ft

NCRP 49– Calculated Requirements for Cross-Table Lateral Wall in Radiographic Room

Using the NCRP 49 attenuation data and recommendations of $W = 1000 \text{ mA-min per wk}$, $U = 1/4$, $T=1$, the new dose limit of $P = 0.02 \text{ mGy (0.002 R) per wk}$, and assuming all exposures are made at 100 kVp , the required barrier thickness is 2.6 mm Pb (1/8 in. or 8 lbs per sq ft).


**Wall Containing Chest Image Receptor
Rad Room**

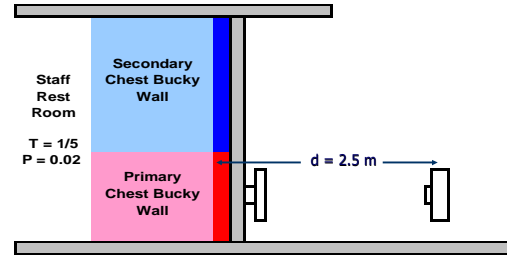

**Wall Containing Chest Image Receptor
Chest Receptor Wall**

- Required thickness $\Leftrightarrow NT/Pd^2$
where:
 - $N = 125 \text{ patients/ week}$
 - $T = 1/5 \text{ (staff rest room)}$
 - $P = 0.02 \text{ mGy/wk}$
 - $d = 2.5 \text{ m}$
- $NT/Pd^2 = 200 \text{ mGy}^{-1} \text{ m}^{-2}$


**Wall Containing Chest Image Receptor
Primary Barrier- Chest Receptor Area**

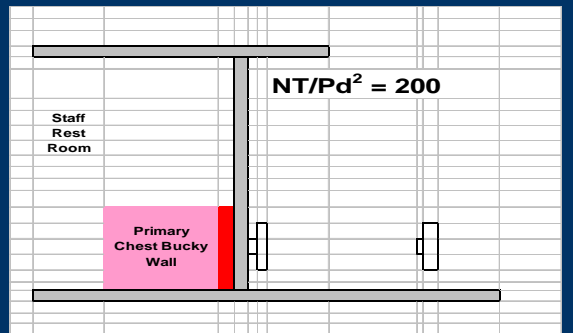

Wall Containing Chest Image Receptor Primary Barrier- Chest Receptor Area

- From Fig 4.5 a, page 54
 - (no pre-shielding)
 - Requires 1.32 mm Pb
- From Fig 4.5 b, page 55
 - (with pre-shielding)
 - Requires 0.50 mm Pb

Wall Containing Chest Image Receptor Secondary Barrier- Chest Receptor Wall

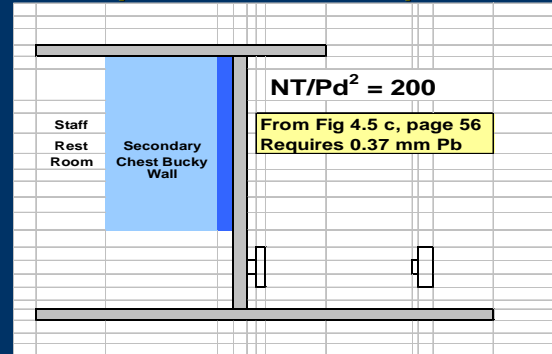


Wall Containing Chest Image Receptor Shielding Required for Entire Wall

Since the primary shielding is greater than the secondary wall requirements, the entire wall can be shielded with the minimum primary requirement.

- No Pre-shielding**
Pb required = 1.32 mm
Specify: 1/16"; 4 lb/sqft
- With Pre-shielding**
Pb required = 0.50 mm
Specify: 1/32"; 2 lb/sqft

Control Wall in the Radiographic Room Secondary Barrier

Controlled Area: $P = 0.1 \text{ mGy/wk}$
 $T = 1$

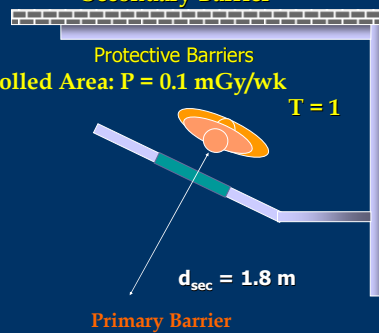

Simplified Graphical Solution
Control Wall in the Radiographic Room

$$NT/Pd^2 = 125 \times 1 / 0.1 \times (1.8)^2 = 386$$


1. Go to page 56, Figure 4.5c "Secondary Wall" curve
2. Look up $NT/Pd^2 = 386$

Pb required = 0.27 mm

Specify: 1/32"; 2 lb/sqft (Minimum)

Recommended Specification: 1/16 inch lead in all control booth walls with lead equivalent windows of at least 1.5 mm.

Floor of the Rad Room

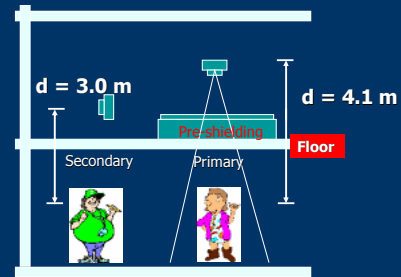

Floor of the Rad Room Primary Barrier Beneath the Rad Table

- Required thickness $\Rightarrow NT/Pd^2$
where:
 - $N = 125$ patients/ week
 - $T = 1$
 - $P = 0.02$ mGy/wk
 - $d = 4.1$ m
- $NT/Pd^2 = 372$ mGy⁻¹ m⁻²

Floor of the Rad Room Primary Barrier Beneath the Rad Table

1. Go to page 58, Fig. 4.6b (Primary, concrete, with pre-shielding)
2. Look up $NT/Pd^2 = 372$

If Specifying: Standard-Weight Concrete:
Minimum Concrete required = 37 mm = 1.5 in.

If Specifying: Light-Weight Concrete:
Minimum Concrete required = 37 mm x 1.3 = 48.1 mm = 1.9 in.

Floor of the Rad Room

Secondary Barrier Calculation for Floor

- Required thickness $\Leftrightarrow NT/Pd^2$
where:
 - $N = 125$ patients/ week
 - $T = 1$
 - $P = 0.02$ mGy/wk
 - $d = 3.0$ m
- $NT/Pd^2 = 694 \text{ mGy}^{-1} \text{ m}^{-2}$

Floor of the Rad Room

Secondary Barrier Calculation for Floor

1. Go to page 59, Fig. 4.6c
(Secondary, concrete)
2. Look up $NT/Pd^2 = 694$

Minimum Concrete required = 33 mm = 1.3 in.

This is less than the 37 mm thickness required for the primary barrier. Thus 37 mm of standard-weight concrete will suffice for the entire floor.

Shielding References

- Simpkin, DJ, Transmission of scatter radiation from computed tomography (CT) scanners determined by a Monte Carlo calculation. *Health Physics* 58(3):363-367, 1990.
- Dixon, RL and Simpkin, DJ. New Concepts for Radiation Shielding of Medical Diagnostic X-ray Facilities. In Proceedings of the 1997 AAPM Summer School.
- NCRP (2005), National Council on Radiation Protection and Measurements. *Structural Shielding Design for Medical X-Ray Imaging Facilities*, NCRP Report #147 (National Council on Radiation Protection and Measurements, Bethesda, Maryland)

Acknowledgement

Multi Slice CT Shielding Slides Courtesy of:

S. Jeff Shephard, M.S., DABR
M.D. Anderson Cancer Center, Houston, TX

Ben Archer, Ph.D, FACR
Baylor College of Medicine, Houston, TX

Multi-Slice Helical CT Shielding

- Larger collimator (slice thickness) settings generate more scatter
 - Offsets advantages of multiple slices per rotation
 - Environmental radiation levels typically increase
- Ceiling and floor deserve close scrutiny

Method

- Calculate the unshielded weekly exposure rate at area of interest
- Find the maximum weekly exposure at 1 m from isocenter and inverse-square this out to the occupied area beyond the barrier.
- Apply traditional barrier thickness calculations to arrive at an answer.
 - Occupancy, permissible dose, attenuation of concrete, etc.

NCRP 147 DLP Method

- Weekly Air Kerma at 1m (K_{sec}^1)

$$K_{sec}^1 \text{ (head)} = K_{head} * DLP$$

$$K_{sec}^1 \text{ (body)} = 1.2 * K_{body} * DLP$$

$$K_{head} = 9 \times 10^{-5} \text{ 1/cm}$$

$$K_{body} = 3 \times 10^{-4} \text{ 1/cm}$$

Use inverse square to find unshielded weekly exposure at barrier from K_{sec}^1

NCRP 147 DLP Method

DLP (Dose-Length Product)

$$= CTDI_{VOL} * L$$

- $CTDI_{VOL} = CTDI_W / \text{Pitch}$

- $CTDI_W = 1/3 \text{ Center } CTDI_{100} + 2/3 \text{ Surface } CTDI_{100} \text{ (mGy)}$

- L = Scan length for average *series* in cm

- Units of mGy-cm

$$= [1/3 CTDI_{100, \text{Center}} + 2/3 CTDI_{100, \text{Surface}}] * L/p$$

NCRP 147 DLP Method

Procedure	CTDI _{Vol} (mGy)	Scan Length (L) (cm)	DLP* (mGy-cm)
Head	60	20	1200
Body	15	35	525
Abdomen	25	25	625
Pelvis	25	20	500
Body (Chest, Abdomen, or Pelvis)			550


* Double the value shown for w/wo contrast


Example - Ceiling Calculation

- 180 Procedures/week
 - 150 Abdomen & Pelvis
 - 30 Head
- 40% w&w/o contrast
- 13' (4.2 m) ceiling height (finished floor to finished floor)
- GE LightSpeed 16

Preliminary Information

- Architectural drawings (Plan view) of exam room, floor above, and floor below
 - Elevation sections through scanner location for floor and ceiling
 - Occupancy factors for floors above and below
 - Two rooms away for possibility that remote areas may be more sensitive than adjacent areas
- Composition of walls, ceilings and floors
 - Materials and thickness
- Scanner placement from vendor
 - Distance from scanner to protected areas beyond barriers

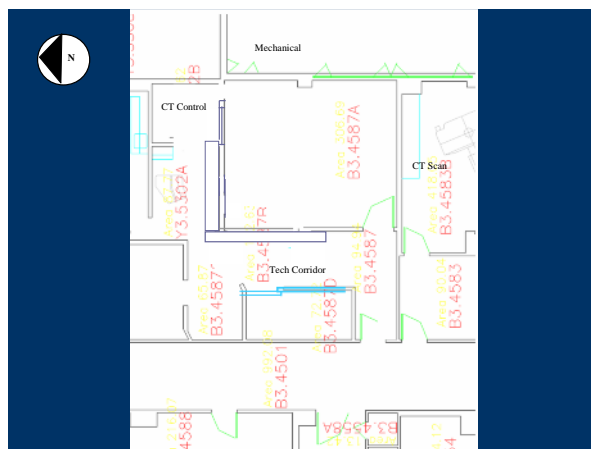

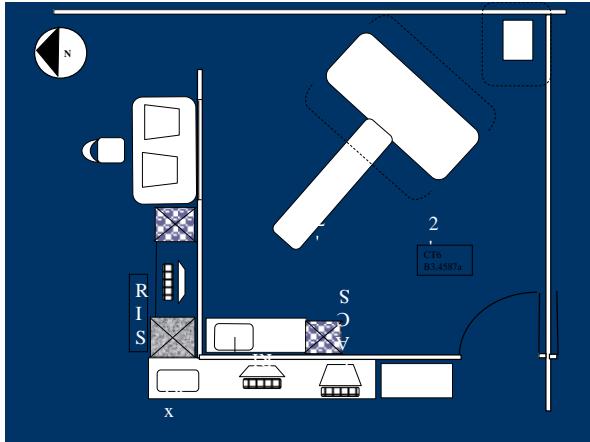

Unshielded Weekly Exposure at Barrier

- Air Kerma/procedure at 1m (K_{sec}^1)
 - 40% w&w/o contrast

$$\begin{aligned}
 K_{sec}^1 \text{ (head)} &= K_{head} * DLP \\
 &= 1.4 * 9 \times 10^{-5} \text{ cm}^{-1} * 1200 \text{ mGy-cm} \\
 &= 4.9 \text{ mGy}
 \end{aligned}$$

$$\begin{aligned}
 K_{sec}^1 \text{ (body)} &= K_{body} * DLP \\
 &= 1.4 * 1.2 * 3 \times 10^{-4} \text{ cm}^{-1} * 550 \\
 &= 41.6 \text{ mGy}
 \end{aligned}$$

Unshielded Weekly Exposure at Barrier

- Weekly Air Kerma (K_{sec}) at Ceiling:
 - 30 head procedures/wk
 - 150 body procedures/wk
 - $D_{sec} = 4.2 \text{ m} + 0.5 \text{ m} - 1 \text{ m} = 3.7 \text{ m}$

$$\begin{aligned}
 K_{sec} \text{ (head)} &= 30 * 4.9 \text{ mGy} * (1\text{m}/3.7\text{m})^2 \\
 &= 0.36 \text{ mGy}
 \end{aligned}$$

$$\begin{aligned}
 K_{sec} \text{ (body)} &= 150 * 41.6 \text{ mGy} * (1\text{m}/3.7\text{m})^2 \\
 &= 3.04 \text{ mGy}
 \end{aligned}$$

Unshielded Weekly Exposure at Barrier

- Weekly Air Kerma (K_{sec}) at Ceiling:

$$K_{sec} \text{ (Total)} = K_{sec} \text{ (head)} + K_{sec} \text{ (body)}$$

$$K_{sec} \text{ (Total)} = 0.36 \text{ mGy} + 3.04 \text{ mGy}$$

$$K_{sec} \text{ (Total)} = 3.40 \text{ mGy}$$

Required Transmission (B)


$$B = \frac{P}{K_{\text{sec}} * T}$$

P = Maximum permissible weekly exposure
T = Occupancy Factor

$$= \frac{0.02 \text{ mGy}}{3.40 \text{ mGy} * 1} = 3.87 \times 10^{-3}$$


Total Shielding Required

Use Simpkin curve fit equations or look up on published attenuation diagrams (NCRP 147 Fig. A-2)

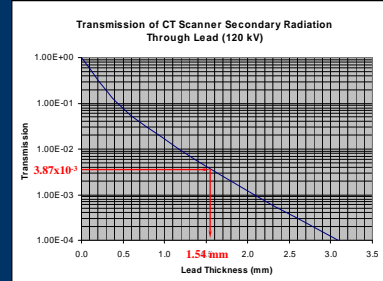




Existing Shielding



- Floors and ceilings
 - Find lead equivalence from documentation of concrete thickness.
 - Find thickness by drilling a test hole and measuring.
 - Always assume light weight concrete, unless proven otherwise (30% less dense than standard density, coefficients used in NCRP 147)

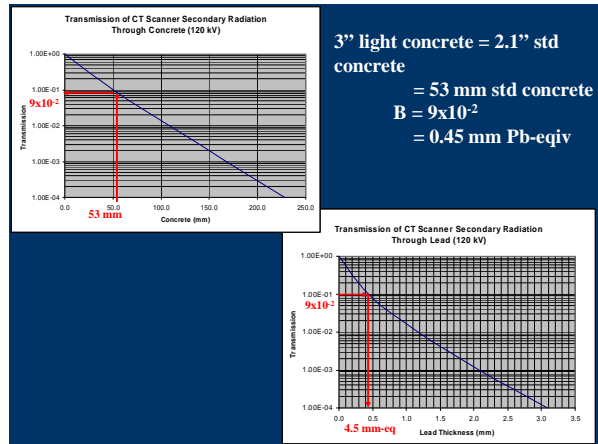

Isodose Map Method

- Assume an isotropic exposure distribution based on the maximum exposure rate in the vendor-supplied exposure distribution plots (approx. 45° to the scanner axis).
- Overestimates shielding needed** in the gantry shadows and the shadows of the patient.

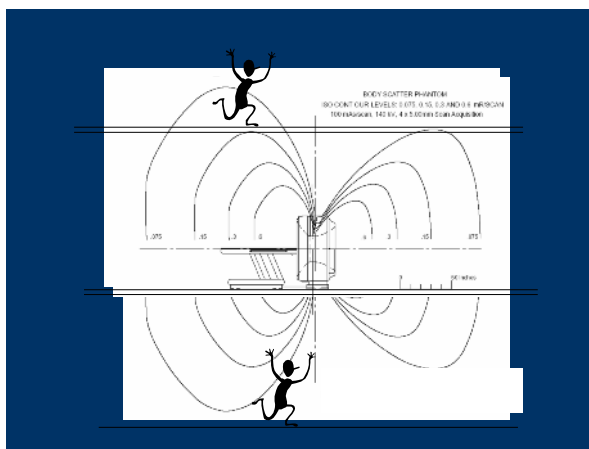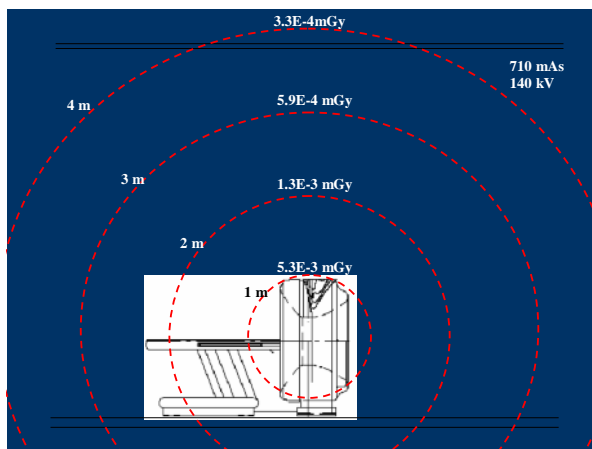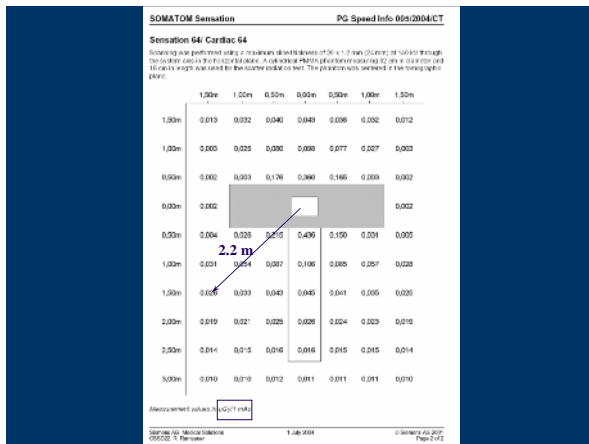

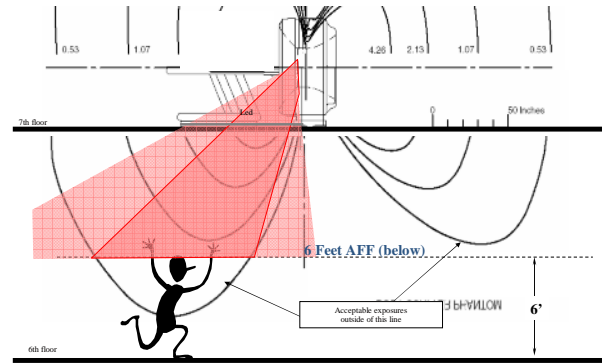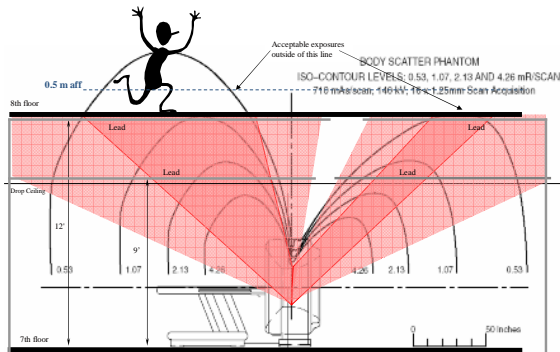

Comparison of Methods						
	DLP		CTDI ₁₀₀		Isodose	
	Head	Body	Head	Body	Head	Body
K ¹ _{sec}	4.9	41.6	0.2	5.0	12	151
Combined Weekly Exposure at Ceiling	3.4 mGy		0.38 mGy		10 mGy	
Add Lead	1/16"		1/32"		3/32"	

Shielding References

- Simpkin, DJ, Transmission of scatter radiation from computed tomography (CT) scanners determined by a Monte Carlo calculation. *Health Physics* 58(3):363-367, 1990.
- Dixon, RL and Simpkin, DJ. New Concepts for Radiation Shielding of Medical Diagnostic X-ray Facilities. In Proceedings of the 1997 AAPM Summer School.
- NCRP (2005), National Council on Radiation Protection and Measurements. *Structural Shielding Design for Medical X-Ray Imaging Facilities*, NCRP Report #147 (National Council on Radiation Protection and Measurements, Bethesda, Maryland)